




LABORATORY LINES INCUBATORS



Kalstein®

A different accompaniment, at your service

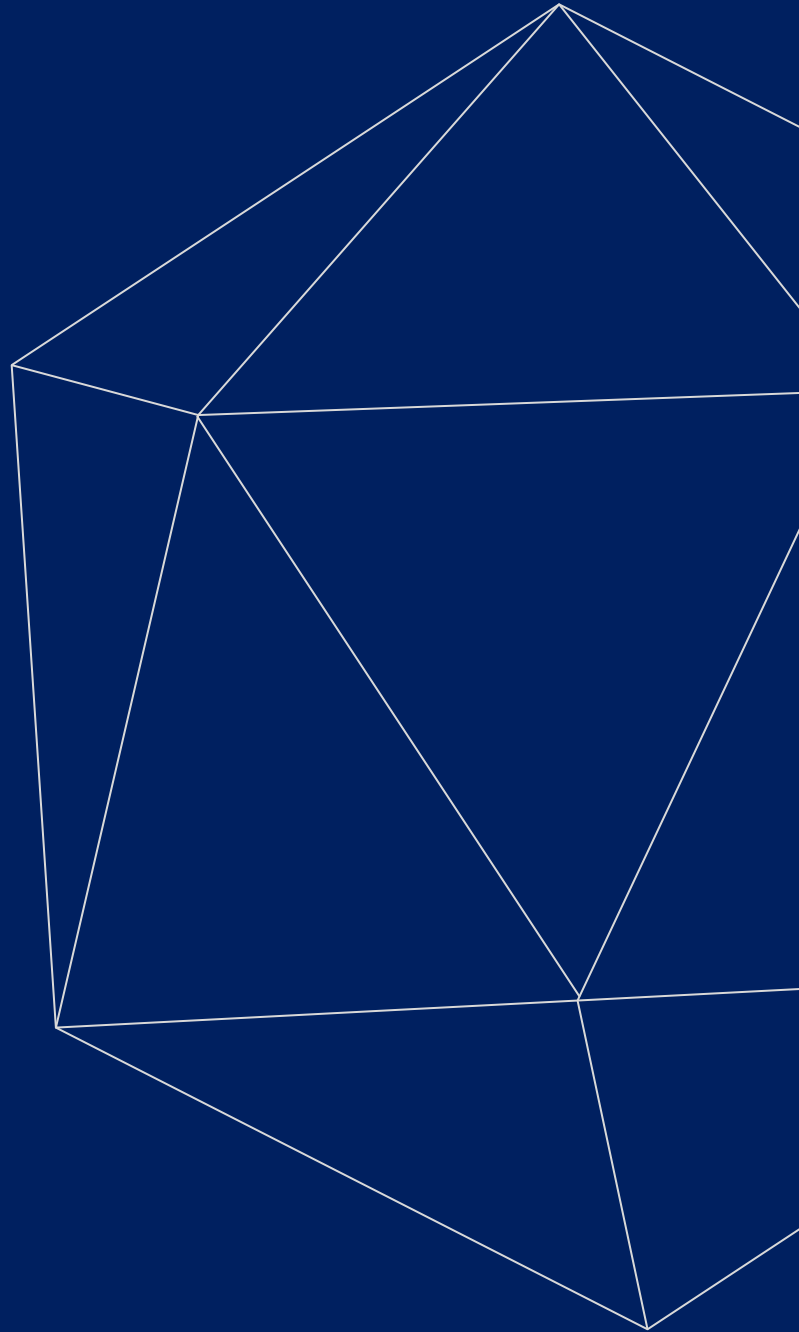

OUR SERVICES

Benefits and Support

In Kalstein France, we take care of the full satisfaction of our customers, that is why we provide value-added services of the highest level based on our experience.



Online Inductions and Trainings

In any part of the world, receive your induction or training from our specialized team of engineers



Quick Response

Our work team is always available to response all your consults or questions, in order to support you in any situation.



#Letsgivemore ♥

Thanks to your purchase, a donation is made to a non-profit foundation that fights against cancer and helps vulnerable communities.



Technical Support

Enjoy of personalized advice for the correct preventive and corrective maintenance of your equipment, thanks to Kalstein's manuals and articles, special catalogues and video tutorials.



Delivery Logistics

We take care of all the necessary logistics for the dispatch of your goods, whether is by sea, land or air.



Kalstein Worldwide

With more than 25 years growing with our customers, Kalstein's multiformat and modern content, is now present in more than 10 countries and increasing.



Biochemical (Mold) Incubator

Series YR02036 (Fluorine free R134A refrigeration)



For environmental protection, medical, drug testing, health and epidemic prevention, livestock, aquaculture and other scientific research and production departments. A special constant temperature equipment for the analysis of the BOD determination of water, bacteria, biochemical, microbiological culture protection, plant cultivation, breeding.

Structure features :

- ✓ The refrigeration system adopts high-quality well-known brand compressors with time-delay protection function and stable operation. It introduces the concept of green environmental protection and uses fluorine-free R134a refrigerant to make the system more energy-saving and efficient and make the environment healthier;
- ✓ Reasonable air duct circulation system to ensure uniform temperature of working room.
- ✓ Mirror stainless steel interior, Shelf spacing and shelf number can be adjusted socket is inbuilt
- ✓ Intelligent microcomputer PID control, touch type parameter setting button, and with timing and over temperature alarm functions
- ✓ Large viewing angle observation window, the inner tank is equipped with lighting device to facilitate observation for YR02036 series)
- ✓ UV light sterilization function, Sterile culture of working area for YR02040 series.

Optional accessories :

- ✓ Programmable controller
- ✓ Test hole
- ✓ Shaker



Model	Inner size (mm) (H*W*D)	Interior Volume (L)	Overall dimension (mm) (H*W*D)	Power (KW)	Temp range (°C)	Temp. Accuracy (°C)	Temp. Fluctuation (°C)
YR02036	500×400×320/	70	1020×540×570	0.35	0~65	±0.5	±1%
YR02037	580×450×380/	100	1100×590×620	0.35	0~65	±0.5	±1%
YR02038	800×500×380/	150	1350×640×660	0.45	0~65	±0.5	±1%
YR02039	1000×550×450/	250	1550×690×660	0.7	0~65	±0.5	±1%
YR02037A	580×450×380/	100	1460*650*700	1	-10~65	±0.5	±1%
YR02038B	800×500×380/	150	1680*700*700	1.2	-10~65	±0.5	±1%
YR02039C	1000×550×450/	250	1880*750*770	1.4	-10~65	±0.5	±1%



#Let'sgive more

Thanks to your purchase, a donation is made to a non-profit foundation that fights against cancer and helps vulnerable communities.



All rights reserved © KALSTEIN France S. A. S.,
Optimum Business Center 450 Rue Baden Powell,
• 34000 Montpellier, France.
Tél: +33 467158849 / +33 680760710/ +33 663810023
<https://kalstein.eu>
KALSTEIN FRANCE, S.A. S



Biochemical (Mold) Incubator

Series YR02036 (Fluorine free R134A refrigeration)



Features

Working principle: electric motor horizontal circulation



#Letsgivemore ♥

Thanks to your purchase, a donation is made to a non-profit foundation that fights against cancer and helps vulnerable communities.




All rights reserved © KALSTEIN France S. A. S.,
Optimum Business Center 450 Rue Baden Powell,
• 34000 Montpellier, France.
Tél: +33 467158849 / +33 680760710/ +33 663810023
<https://kalstein.eu>
KALSTEIN FRANCE, S.A. S



Biochemical (Mold) Incubator Series YR02036 (Fluorine free R134A refrigeration)



⊕ Stainless steel liner

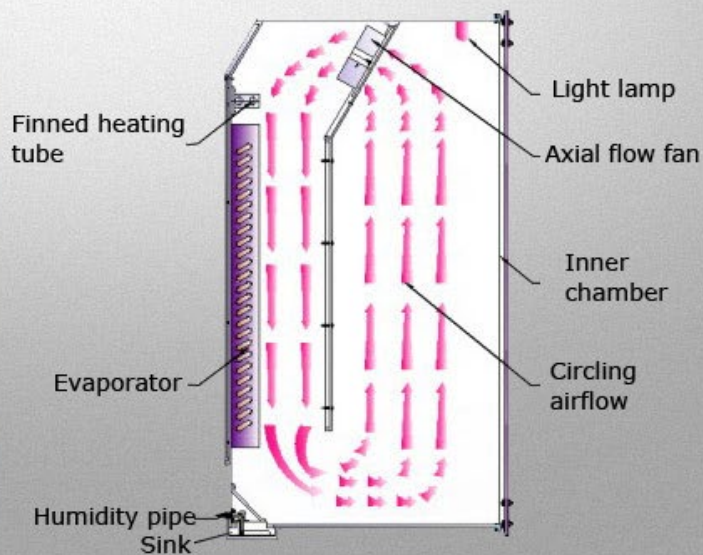

⊕ ALLFLOW™ perfect Air Current Cycling



⊕ Without fluorine R134a refrigeration compressor



⊕ Lighting system and BOD socket



#Let'sgive**more** ♥

Thanks to your purchase, a donation is made to a non-profit foundation that fights against cancer and helps vulnerable communities.



All rights reserved © KALSTEIN France S. A. S.,
Optimum Business Center 450 Rue Baden Powell,
• 34000 Montpellier, France.
Tél: +33 467158849 / +33 680760710/ +33 663810023
<https://kalstein.eu>
KALSTEIN FRANCE, S.A. S



Biochemical (Mold) Incubator



Series YR02036 (Fluorine free R134A refrigeration)



⊕ Intelligent digital display

Operation chart:



1

Press the power switch
(Specific operation can refer to the description of the instructions)






2

Press the set button 
Press once time to set the temperature, press two time to set the time, press three times to confirm and exit



3

Set time and temperature
( Select shift key,  add and subtract key,  lighting key)



4

Set up to complete
Automatic operation of the machine, arrive time will be automatic alarm, and then stop running.

Optional accessories:



1. Test hole



2. The governor multi-use for oscillator



3. Intelligent programmable temperature controller



 **Kalstein**[®]
A different accompaniment, at your service

#LETSGIVEMORE 

**WITH THE
ACQUISITION
OF A KALSTEIN
EQUIPMENT**

YOU MAKE A CONTRIBUTION TO:



Fundación
Jacinto Convit

<https://www.jacintoconvit.org/>



<https://onetreepanted.org>



Humatem

Solidarité et coopération
biomédicale


<https://www.humatem.org/>



**MANIAPURE
FOUNDATION USA**

<https://www.maniapurefoundation.org/>

Articles Kalstein
(Social Networks and Official Websites)

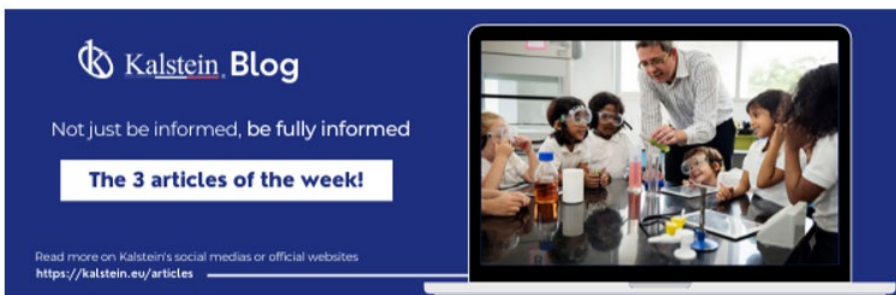

Kalstein Blog
Not just be informed, be fully informed
The 3 articles of the week!
Read more on Kalstein's social medias or official websites
<https://kalstein.eu/articles>

1
Read this ARTICLE NOW



How does the climate affect your microscope?
A microscope is an instrument that allows observing and measuring objects that are too small to be seen with the naked eye. This is achieved through a complex system of lenses and mirrors, which form

Kalstein Blog
Article of the week:
How does the climate affect your microscope?
Read more on Kalstein's social medias or official websites
<https://kalstein.eu/articles>

2
Read this ARTICLE NOW



What is the importance of the laboratory refrigerator?
A refrigerator in a laboratory is one of the most important pieces of equipment. It is used to maintain, in a controlled environment (refrigerated space), various fluids and substances.

Kalstein Blog
Article of the week:
What is the importance of the laboratory refrigerator?
Read more on Kalstein's social medias or official websites
<https://kalstein.eu/articles>

Kalstein Training Videos
(Youtube Channel and Official Websites)



Kalstein France - Jacinto Convit Foundation / Alliances
23 views • 3 years ago



Kalstein's Microscopes
445 views • 1 week ago



MONTPELLIER - FRANCE
HEADQUARTER
Optimum Business Center,
450 Rue Baden Powell,
34000 Montpellier - France.

Fax: +33 (0) 4 67 15 88 49
Tif: +33 (0) 6 80 76 07 10
sales@kalstein.eu

<https://kalstein.net/>
<https://www.kalstein.fr/>

Official Websites for Latin America.

