




LABORATORY SECTOR

BIOSAFETY CABINETS



Kalstein®

A different accompaniment, at your service

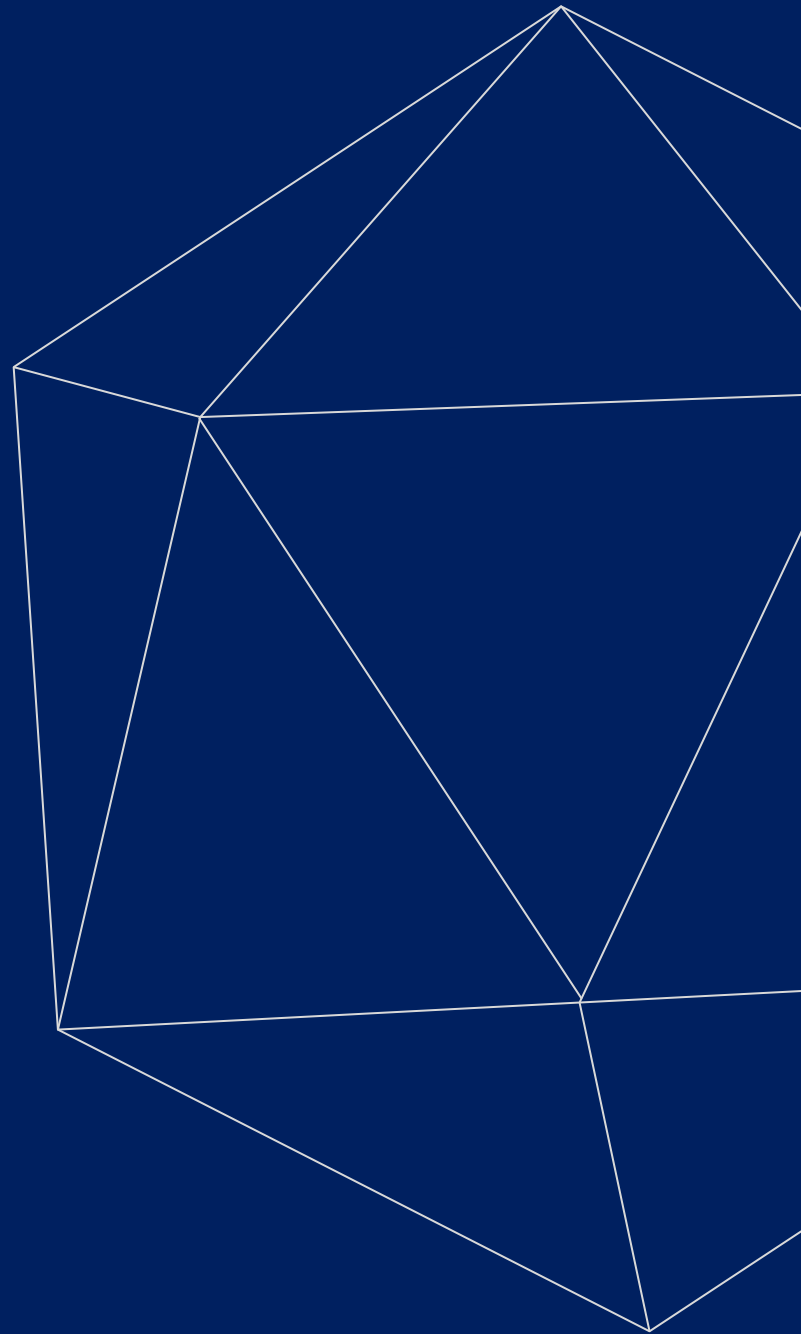



Kalstein®

A different accompaniment, at your service

OUR SERVICES

Benefits and Support

In Kalstein France, we take care of the full satisfaction of our customers, that is why we provide value-added services of the highest level based on our experience.



Online Inductions and Trainings

In any part of the world, receive your induction or training from our specialized team of engineers



Quick Response

Our work team is always available to response all your consults or questions, in order to support you in any situation.



#Letsgivemore

Thanks to your purchase, a donation will be made to a non-profit foundation that fights against breast cancer and helps most vulnerable communities.



Technical Support

Enjoy of personalized advice for the correct preventive and corrective maintenance of your equipment, thanks to Kalstein's manuals and articles, special catalogues and video tutorials.



Delivery Logistics

We take care of all the necessary logistics for the dispatch of your goods, whether is by sea, land or air.



Kalstein Worldwide

With more than 25 years growing with our customers, Kalstein's multiformat and modern content, is now present in more than 10 countries and increasing.



KALSTEIN® Biosafety Cabinets
Biological Safety Cabinet
Model YR05708 // YR05710



#Letsgivemore ♥

Thanks to your purchase, a donation will be made to a non-profit foundation that fights against breast cancer and helps most vulnerable communities.



All rights reserved © KALSTEIN France S. A. S.,
2 Rue Jean Lantier • 75001 Paris •
+33 1 78 95 87 89 / +33 6 80 76 07 10 •
<https://kalstein.eu>
KALSTEIN FRANCE, S. A. S



KALSTEIN® Biosafety Cabinets Biological Safety Cabinet Model YR05708 // YR05710



Product Description:

YR05708 // YR05710 biological safety cabinet Class II biological safety cabinet/biological safety cabinet is necessary in animalcule lab, especially in the condition needing protection measure, such as medicine, pharmacy, scientific research and so on.

Features:

- ✓ Air curtain isolation design prevents internal and external cross-contamination, 100% air flow is discharged outside, negative pressure vertical laminar flow, no need to install pipes
- ✓ The glass door can be moved up and down, can be positioned arbitrarily, is easy to operate, and can be completely closed for sterilization, and the positioning height limit alarm prompts.
- ✓ The power output socket in the work area is equipped with a waterproof socket and a sewage interface to provide great convenience for the operator
- ✓ A special filter is installed at the exhaust air to control emission pollution
- ✓ The working environment is made of high-quality 304 stainless steel, which is smooth, seamless, and has no dead ends. It can be easily and thoroughly disinfected and can prevent corrosion by corrosive agents and disinfectants
- ✓ It adopts LED LCD panel control and built-in UV lamp protection device, which can only be opened when the safety door is closed
- ✓ With DOP detection port, built-in differential pressure gauge.
- ✓ 10° tilt angle, in line with the human body design concept



#Letsgivemore ♥

Thanks to your purchase, a donation will be made to a non-profit foundation that fights against breast cancer and helps most vulnerable communities.



All rights reserved © KALSTEIN France S. A. S.,
2 Rue Jean Lantier • 75001 Paris •
+33 1 78 95 87 89 / +33 6 80 76 07 10 •
<https://kalstein.eu>
KALSTEIN FRANCE, S. A. S



KALSTEIN® Biosafety Cabinets Biological Safety Cabinet Model YR05708 // YR05710



Technical Specifications:

Model	YR05708	YR05709	YR05710
Cleanliness	Class 100@≥0.5µm (US Federal 209E)		
Number of colonies	≤0.5pcs/dish·hour (Φ90mm culture plate)		
average wind speed	Inside the door	0.38±0.025m/s	
	middle	0.26±0.025m/s	
	Inside	0.27±0.025m/s	
Front suction wind speed	0.55m±0.025m/s		
noise	≤65dB(A)		
Voltage	AC single-phase 220V/50Hz		
Maximum power consumption	500W	600W	700W
Weight	210KG	250KG	270KG
Working area W₁×D₁×H₁	1040×650×620	1340×650×620	1640×650×620
Overall dimension W×D×H	1200×800×2270	1500×800×2270	1800×800×2270
High efficiency filter specification and quantity	980×490×50×1	1280×490×50×1	1580×490×50×1
	520×380×70×1	820×380×70×1	1120×380×70×1
Specification and quantity of fluorescent lamp/ultraviolet lamp	20W×2/20W×1	30W×2/30W×1	30W×2/30W×1



#Letsgivemore

Thanks to your purchase, a donation will be made to a non-profit foundation that fights against breast cancer and helps most vulnerable communities.



All rights reserved © KALSTEIN France S. A. S.,
2 Rue Jean Lantier • 75001 Paris •
+33 1 78 95 87 89 / +33 6 80 76 07 10 •
<https://kalstein.eu>
KALSTEIN FRANCE, S. A. S





 **Kalstein**[®]
A different accompaniment, at your service

#LETSGIVEMORE 

**WITH THE
ACQUISITION
OF A KALSTEIN
EQUIPMENT**

YOU MAKE A CONTRIBUTION TO:



Fundación
Jacinto Convit

<https://www.jacintoconvit.org/>



<https://onetreepanted.org>



Humatem

Solidarité et coopération
biomédicale

<https://www.humatem.org/>



**MANIAPURE
FOUNDATION USA**

<https://www.maniapurefoundation.org/>

Articles Kalstein
(Social Networks and Official Websites)



Kalstein Blog
Not just be informed, be fully informed

The 3 articles of the week!

Read more on Kalstein's social medias or official websites
<https://kalstein.eu/articles>

1
Read this ARTICLE NOW



How does the climate affect your microscope?

microscope is an instrument that allows observing and measuring small objects. This is possible through the use of lenses or a system of mirrors that collect and focus light from the object being observed.

Kalstein Blog

Article of the week:
How does the climate affect your microscope?

Read more on Kalstein's social medias or official websites
<https://kalstein.eu/articles>

2
Read this ARTICLE NOW



What is the importance of the laboratory refrigerator?

Refrigerator is a device that is one of the most common equipment that is used to maintain, in a controlled environment (refrigerated space), various fluids and substances.

Kalstein Blog

Article of the week:
What is the importance of the laboratory refrigerator?

Read more on Kalstein's social medias or official websites
<https://kalstein.eu/articles>

Kalstein Training Videos
(Youtube Channel and Official Websites)



Kalstein France - Jacinto Convit Foundation / Allianc...
23 views • 3 years ago



Kalstein's Microscopes
445 views • 1 week ago



PARIS - FRANCE
HEADQUARTER
2 Rue Jean Lantier
75001 Paris - France

Tlf: +33 (0) 01 78 95 87 89
Cel: +33 (0) 6 80 76 07 10
sales@kalstein.eu

<https://kalstein.eu/>
<https://www.kalstein.fr/>

