



LABORATORY SECTOR

REFRIGERATORS AND FREEZERS



Kalstein®

A different accompaniment, at your service





Kalstein®

A different accompaniment, at your service

OUR SERVICES

Benefits and Support

In Kalstein France, we take care of the full satisfaction of our customers, that is why we provide value-added services of the highest level based on our experience.



Online Inductions and Trainings

In any part of the world, receive your induction or training from our specialized team of engineers



Quick Response

Our work team is always available to response all your consults or questions, in order to support you in any situation.



#Letsgivemore

Thanks to your purchase, a donation will be made to a non-profit foundation that fights against breast cancer and helps most vulnerable communities.



Technical Support

Enjoy of personalized advice for the correct preventive and corrective maintenance of your equipment, thanks to Kalstein's manuals and articles, special catalogues and video tutorials.



Delivery Logistics

We take care of all the necessary logistics for the dispatch of your goods, whether is by sea, land or air.



Kalstein Worldwide

With more than 25 years growing with our customers, Kalstein's multiformat and modern content, is now present in more than 10 countries and increasing.



KALSTEIN® Refrigerators and Freezers
Blood Bank Freezer With 12 Pcs 304 Stainless Drawers

Model YR05292



#Letsgivemore ♥

Thanks to your purchase, a donation will be made to a non-profit foundation that fights against breast cancer and helps most vulnerable communities.



All rights reserved © KALSTEIN France S. A. S.,
2 Rue Jean Lantier • 75001 Paris •
+33 1 78 95 87 89 / +33 6 80 76 07 10 •
<https://kalstein.eu>
KALSTEIN FRANCE, S. A. S



KALSTEIN® Refrigerators and Freezers

Blood Bank Freezer With 12 Pcs 304 Stainless Drawers

Model YR05292



Product Description:

1008L blood bank refrigerator with 12 pcs 304 stainless steel drawers , temperature $4\pm 1^{\circ}\text{C}$.

PROMED Blood Bank Refrigerator are Specially Designed to Store Blood Store Blood Bags and Blood Derived Products:

For 450 ml blood bags storage: 576 units

Features:

- ✓ Accurate microprocessor controller to ensure consistent cabinet temperature $4\pm 1^{\circ}\text{C}$.
- ✓ Specially designed forced-air cooling system to ensure frost-free cabinet and good temperature uniformity.
- ✓ Three layer toughened glass door with heater to avoid condensation.
- ✓ 304 stainless steel inner chamber for friendly using and easy cleaning.
- ✓ Wide available voltage range: 187V~242V

Safety:

- ✓ Safety door lock and controller setting password to prevent unauthorized access.
- ✓ Build-in back-up battery to display temperature and to operate audible and visual alarm system for up to 8 hours without AC power.
- ✓ Inner doors and auto on/off cooling fan to prevent cold air from leaking out when the door is opening.
- ✓ Visual and audible alarm system: High/low temperature, power failure, sensor error, door ajar, backup battery malfunction
- ✓ Power failure protection: turn on delay of the cooling system after power failure.
- ✓ Controller error protection: the cooling system will remember the normal working cycle and keep working depending on this memory when controller or two sensor failed.
- ✓ Temperature printer: print current inside temperature automatically per 15 minutes.

Humanized Design:

- ✓ 4 casters with stopper for easy moving and fixation.
- ✓ 12 stainless steel drawers for easy using.
- ✓ Special designed service-free drainage system ensure no water in the cabinet.
- ✓ Standard 25mm test hole for easy connecting to temperature monitoring systems.
- ✓ Digital display with big screen, easy to check



#Letsgivemore 

Thanks to your purchase, a donation will be made to a non-profit foundation that fights against breast cancer and helps most vulnerable communities.



All rights reserved © KALSTEIN France S. A. S.,
2 Rue Jean Lantier • 75001 Paris •
+33 1 78 95 87 89 / +33 6 80 76 07 10 •
<https://kalstein.eu>
KALSTEIN FRANCE, S. A. S



KALSTEIN® Refrigerators and Freezers

Blood Bank Freezer With 12 Pcs 304 Stainless Drawers

Model YR05292



Technical Description:

Model	YR05292
Performance	
Temperature Range	4°C±1°C
Ambient Temperature	10~32°C
Controller	Microprocessor
Display	Digital Display
Refrigeration System	
Compressor	1
Refrigerant	R134a
Cooling Type	Forced-air cooling
Defrost Type	Frost free
Alarm	
High/Low Temperature	Y
Sensor Error	Y
Power Failure	Y
Door ajar	Y



#Letsgivemore

Thanks to your purchase, a donation will be made to a non-profit foundation that fights against breast cancer and helps most vulnerable communities.



All rights reserved © KALSTEIN France S. A. S.,
2 Rue Jean Lantier • 75001 Paris •
+33 1 78 95 87 89 / +33 6 80 76 07 10 •
<https://kalstein.eu>
KALSTEIN FRANCE, S. A. S



KALSTEIN® Refrigerators and Freezers

Blood Bank Freezer With 12 Pcs 304 Stainless Drawers

Model YR05292



Technical Description:

Accessories	
Caster	4
Lamp	LED
Drawer	12
Test hole	1/25mm
Temperature printer	1
Chart recorder	Optional
Electrical Data	
Power Supply(V/Hz)	220/50~60; 115/60
Rated Power(W)	513
Electrical Current(A)	2.74
Power Consumption	8.53 kWh/24h
Dimensions	
Capacity	1008L
Interior Size(W*D*H)	1100*684*1325mm
Exterior Size(W*D*H)	1240*860*1895mm
Package Size(W*D*H)	1300*960*2130mm
N.t./G.t.(approx)	265/335kg
Loading Qty(20'/40'/40'H)	8/16/16



#Letsgivemore

Thanks to your purchase, a donation will be made to a non-profit foundation that fights against breast cancer and helps most vulnerable communities.



All rights reserved © KALSTEIN France S. A. S.,
2 Rue Jean Lantier • 75001 Paris •
+33 1 78 95 87 89 / +33 6 80 76 07 10 •
<https://kalstein.eu>
KALSTEIN FRANCE, S. A. S





 **Kalstein**[®]
A different accompaniment, at your service

#LETSGIVEMORE 

**WITH THE
ACQUISITION
OF A KALSTEIN
EQUIPMENT**

YOU MAKE A CONTRIBUTION TO:



**Fundación
Jacinto Convit**

<https://www.jacintoconvit.org/>



<https://onetreepanted.org>



Humatem

Solidarité et coopération
biomédicale


<https://www.humatem.org/>



**MANIAPURE
FOUNDATION USA**

<https://www.maniapurefoundation.org/>

Articles Kalstein
(Social Networks and Official Websites)

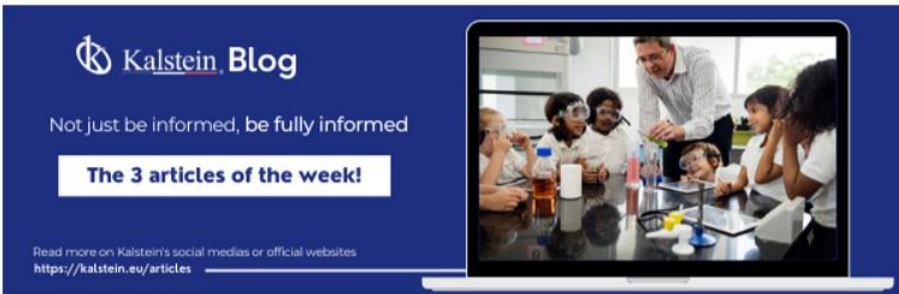

Kalstein Blog
Not just be informed, be fully informed

The 3 articles of the week!

Read more on Kalstein's social medias or official websites
<https://kalstein.eu/articles>



1
Read this ARTICLE NOW

How does the climate affect your microscope?

Kalstein Blog

Article of the week:
How does the climate affect your microscope?

Read more on Kalstein's social medias or official websites
<https://kalstein.eu/articles>



2
Read this ARTICLE NOW

What is the importance of the laboratory refrigerator?

Kalstein Blog

Article of the week:
What is the importance of the laboratory refrigerator?

Read more on Kalstein's social medias or official websites
<https://kalstein.eu/articles>

Kalstein Training Videos
(Youtube Channel and Official Websites)



Kalstein France - Jacinto Convit Foundation / Allianc...
23 views • 3 years ago



Kalstein's Microscopes
445 views • 1 week ago



PARIS - FRANCE
HEADQUARTER
2 Rue Jean Lantier
75001 Paris - France

Tlf: +33 (0) 01 78 95 87 89
Cel: +33 (0) 6 80 76 07 10
sales@kalstein.eu

<https://kalstein.eu/>
<https://www.kalstein.fr/>

