




# LABORATORY SECTOR

## HOTPLATES AND HOTPLATE STIRRERS



**Kalstein**®

A different accompaniment, at your service

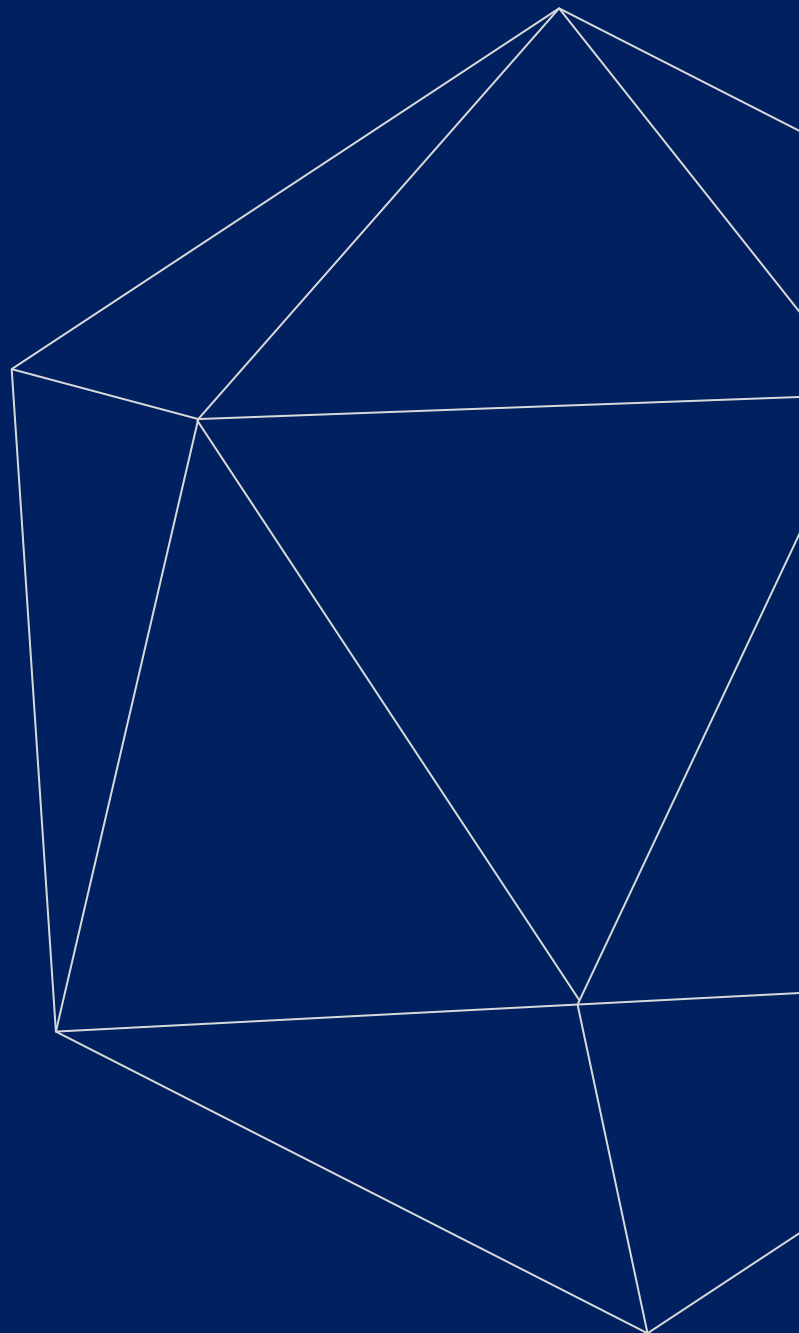



# Kalstein®

A different accompaniment, at your service

OUR SERVICES

## Benefits and Support

In Kalstein France, we take care of the full satisfaction of our customers, that is why we provide value-added services of the highest level based on our experience.



### Online Inductions and Trainings

In any part of the world, receive your induction or training from our specialized team of engineers



### Quick Response

Our work team is always available to response all your consults or questions, in order to support you in any situation.



### #Letsgivemore ♥

Thanks to your purchase, a donation will be made to a non-profit foundation that fights against breast cancer and helps most vulnerable communities.



### Technical Support

Enjoy of personalized advice for the correct preventive and corrective maintenance of your equipment, thanks to Kalstein's manuals and articles, special catalogues and video tutorials.



### Delivery Logistics

We take care of all the necessary logistics for the dispatch of your goods, whether is by sea, land or air.



### Kalstein Worldwide

With more than 25 years growing with our customers, Kalstein's multiformat and modern content, is now present in more than 10 countries and increasing.



**KALSTEIN® Hotplates and Hotplate Stirrers**  
**Industrial Hotplate Magnetic Stirrer Bar With Heating**  
**Model YR02939**



**#Letsgivingore** ♥

Thanks to your purchase, a donation will be made to a non-profit foundation that fights against breast cancer and helps most vulnerable communities.



All rights reserved © KALSTEIN France S. A. S.,  
2 Rue Jean Lantier • 75001 Paris •  
+33 1 78 95 87 89 / +33 6 80 76 07 10 •  
<https://kalstein.eu>  
**KALSTEIN FRANCE, S. A. S**



# KALSTEIN® Hotplates and Hotplate Stirrers

## Industrial Hotplate Magnetic Stirrer Bar With Heating

### Model YR02939



#### Product Description:

Magnetic stirrers are often used in chemistry and biology. They are preferred over gear-driven motorized stirrers because they are quieter, more efficient, and have no moving external parts to break or wear out (other than the simple bar magnet itself). Because of its small size, a stirring bar is more easily cleaned and sterilized than other stirring devices. They do not require lubricants which could contaminate the reaction vessel and the product. They can be used inside hermetically closed vessels or systems, without the need for complicated rotary seals. Magnetic stirrers may also include a hot plate or some other means for heating the liquid.

#### Product Features:

- ✓ Very compact and simple to operate
- ✓ Digital Display temperature and speed, Heating temperature: Rt-380 Degree and mixing speed :100-1800 RPM
- ✓ Heat and corrosion resistant aluminum alloy housing for durability
- ✓ Ceramic coated aluminum heating plate for great heat distribution and corrosion resistance
- ✓ Specially designed housing prevents spill from damaging the electronic components inside
- ✓ Secondary heat insulation plate prevents heat loss and eliminates heat damage to electronic components
- ✓ Adjustable temperature and rotation speed



#Letsgivemore ♥

Thanks to your purchase, a donation will be made to a non-profit foundation that fights against breast cancer and helps most vulnerable communities.



All rights reserved © KALSTEIN France S. A. S.,  
2 Rue Jean Lantier • 75001 Paris •  
+33 1 78 95 87 89 / +33 6 80 76 07 10 •  
<https://kalstein.eu>  
KALSTEIN FRANCE, S. A. S.



# KALSTEIN® Hotplates and Hotplate Stirrers

## Industrial Hotplate Magnetic Stirrer Bar With Heating

### Model YR02939



#### Technical Description:

|  |  |
|--|--|
| <b>Model</b>                           | YR02939                                  |
| <b>Volume</b>                          | 10L                                      |
| <b>Heating Plate (mm)</b>              | 230*230                                  |
| <b>Heating Power(W)</b>                | 750                                      |
| <b>Heating Plate Material</b>          | Ceramic Coated Aluminum (White or Black) |
| <b>Temperature Range</b>               | Rt-380°C                                 |
| <b>Mixing Speed</b>                    | 100-1800RPM                              |
| <b>Temperature Control</b>             | PID Auto Control                         |
| <b>Temperature Display (°C )</b>       | Digital LCD Display                      |
| <b>Temperature Accuracy</b>            |  |
| <b>Thermocouple Probe(Sensor)</b>      | PT100, with Inner and Outer Thermocouple |
| <b>Permissible Working Temperature</b> | Rt to 40°C                               |
| <b>Voltage</b>                         | 220V/ 110V Optional                      |
| <b>Dimensions (mm)</b>                 | 500 x 400 x 350                          |
| <b>Weight (Kg)</b>                     | 7.5                                      |



**#Letsgivemore** ♥

Thanks to your purchase, a donation will be made to a non-profit foundation that fights against breast cancer and helps most vulnerable communities.



All rights reserved © KALSTEIN France S. A. S.,  
2 Rue Jean Lantier • 75001 Paris •  
+33 1 78 95 87 89 / +33 6 80 76 07 10 •  
<https://kalstein.eu>  
KALSTEIN FRANCE, S. A. S





 **Kalstein**<sup>®</sup>  
A different accompaniment, at your service

#LETSGIVEMORE 

**WITH THE  
ACQUISITION  
OF A KALSTEIN  
EQUIPMENT**

YOU MAKE A CONTRIBUTION TO:



**Fundación  
Jacinto Convit**

<https://www.jacintoconvit.org/>



<https://onetreepanted.org>



**Humatem**

Solidarité et coopération  
biomédicale


<https://www.humatem.org/>



**MANIAPURE  
FOUNDATION USA**

<https://www.maniapurefoundation.org/>

**Articles Kalstein**  
(Social Networks and Official Websites)

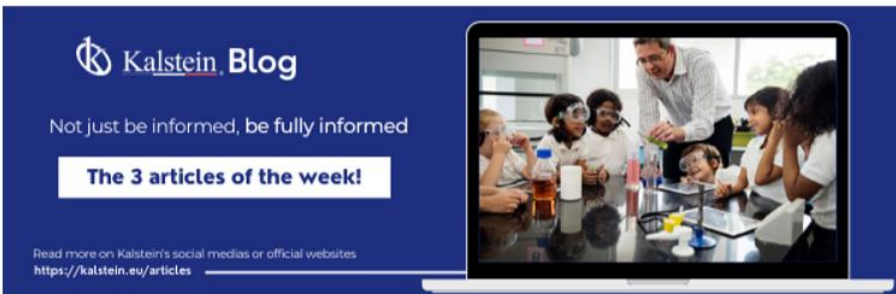

**Kalstein Blog**  
Not just be informed, be fully informed  
**The 3 articles of the week!**  
Read more on Kalstein's social medias or official websites  
<https://kalstein.eu/articles>



**1**  
Read this ARTICLE NOW  
**How does the climate affect your microscope?**  
microscope in a laboratory. One of the most important pieces of equipment that needs to maintain in a controlled environment (refrigerated space), without fluctuations in temperature.  
**Kalstein Blog**  
**Article of the week:**  
How does the climate affect your microscope?  
Read more on Kalstein's social medias or official websites  
<https://kalstein.eu/articles>



**2**  
Read this ARTICLE NOW  
**What is the importance of the laboratory refrigerator?**  
Refrigerator in a laboratory is one of the most important pieces of equipment that needs to maintain in a controlled environment (refrigerated space), without fluctuations in temperature.  
**Kalstein Blog**  
**Article of the week:**  
What is the importance of the laboratory refrigerator?  
Read more on Kalstein's social medias or official websites  
<https://kalstein.eu/articles>

**Kalstein Training Videos**  
(Youtube Channel and Official Websites)



**Kalstein France - Jacinto Convit Foundation / Allianc...**  
23 views • 3 years ago



**Kalstein's Microscopes**  
445 views • 1 week ago



PARIS - FRANCE  
**HEADQUARTER**  
2 Rue Jean Lantier  
75001 Paris - France

**Tlf:** +33 (0) 01 78 95 87 89  
**Cel:** +33 (0) 6 80 76 07 10  
[sales@kalstein.eu](mailto:sales@kalstein.eu)

<https://kalstein.eu/>  
<https://www.kalstein.fr/>

