



LABORATORY SECTOR

HOTPLATES AND HOTPLATE STIRRERS



Kalstein®

A different accompaniment, at your service





Kalstein®

A different accompaniment, at your service

OUR SERVICES

Benefits and Support

In Kalstein France, we take care of the full satisfaction of our customers, that is why we provide value-added services of the highest level based on our experience.



Online Inductions and Trainings

In any part of the world, receive your induction or training from our specialized team of engineers



Quick Response

Our work team is always available to response all your consults or questions, in order to support you in any situation.



#Letsgivemore

Thanks to your purchase, a donation will be made to a non-profit foundation that fights against breast cancer and helps most vulnerable communities.



Technical Support

Enjoy of personalized advice for the correct preventive and corrective maintenance of your equipment, thanks to Kalstein's manuals and articles, special catalogues and video tutorials.



Delivery Logistics

We take care of all the necessary logistics for the dispatch of your goods, whether is by sea, land or air.




Kalstein Worldwide

With more than 25 years growing with our customers, Kalstein's multiformat and modern content, is now present in more than 10 countries and increasing.



KALSTEIN® Hotplates and Hotplate Stirrers Magnetic Stirrer With Hotplate Model YR02951 // YR02952

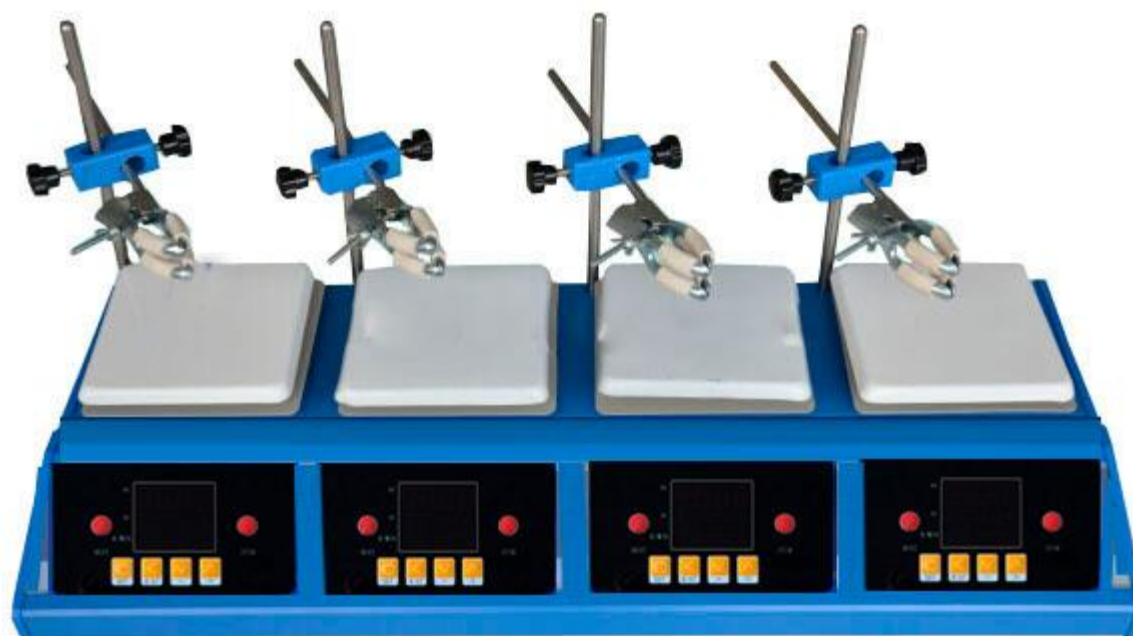


#Letsgivingmore ♥

Thanks to your purchase, a donation will be made to a non-profit foundation that fights against breast cancer and helps most vulnerable communities.



All rights reserved © KALSTEIN France S. A. S.,
2 Rue Jean Lantier • 75001 Paris •
+33 1 78 95 87 89 / +33 6 80 76 07 10 •
<https://kalstein.eu>
KALSTEIN FRANCE, S. A. S

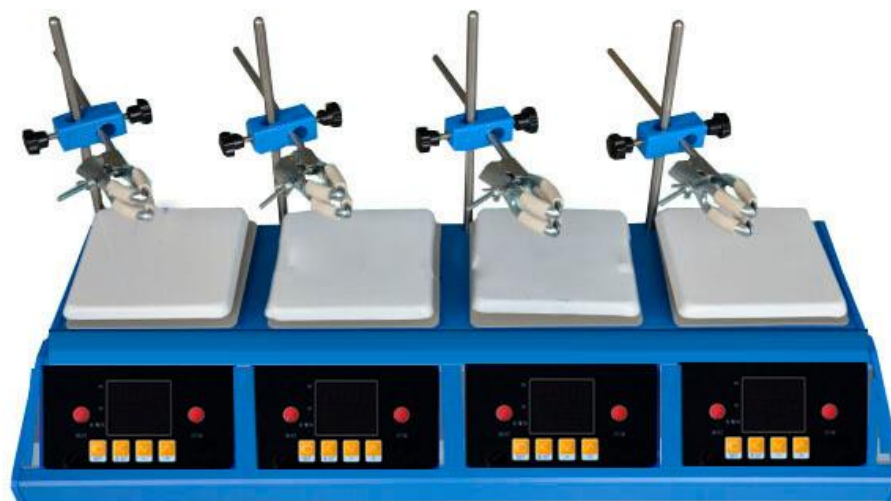

KALSTEIN® Hotplates and Hotplate Stirrers**Magnetic Stirrer With Hotplate****Model YR02951 // YR02952****Product Description:**

Magnetic stirrers are often used in chemistry and biology. They are preferred over gear-driven motorized stirrers because they are quieter, more efficient, and have no moving external parts to break or wear out (other than the simple bar magnet itself). Because of its small size, a stirring bar is more easily cleaned and sterilized than other stirring devices. They do not require lubricants which could contaminate the reaction vessel and the product. They can be used inside hermetically closed vessels or systems, without the need for complicated rotary seals. Magnetic stirrers may also include a hot plate or some other means for heating the liquid.

Product Features:

- ✓ Very compact and simple to operate
- ✓ Digital Display temperature and speed, Heating temperature: Rt-380 Degree and mixing speed :100-1800 RPM
- ✓ Heat and corrosion resistant aluminum alloy housing for durability
- ✓ Ceramic coated aluminum heating plate for great heat distribution and corrosion resistance
- ✓ Specially designed housing prevents spill from damaging the electronic components inside
- ✓ Secondary heat insulation plate prevents heat loss and eliminates heat damage to electronic components
- ✓ Adjustable temperature and rotation speed

**#Letsgivemore** ♥

Thanks to your purchase, a donation will be made to a non-profit foundation that fights against breast cancer and helps most vulnerable communities.




All rights reserved © KALSTEIN France S. A. S.,
2 Rue Jean Lantier • 75001 Paris •
+33 1 78 95 87 89 / +33 6 80 76 07 10 •
<https://kalstein.eu>
KALSTEIN FRANCE, S. A. S



KALSTEIN® Hotplates and Hotplate Stirrers

Magnetic Stirrer With Hotplate

Model YR02951 // YR02952

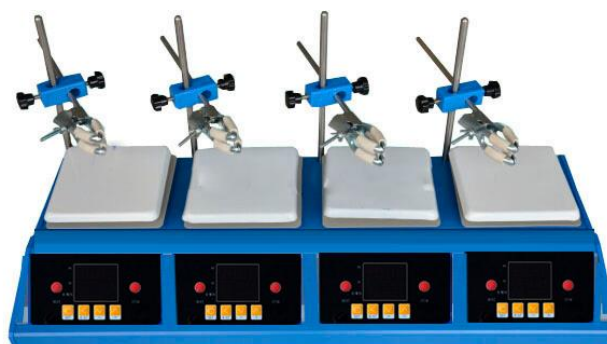

Technical Description:

Model	YR02951	YR02952
Upper plate material	Ceramic White	
Remove from	0 ~ 2500rpm	
Remove speed accuracy	± 1 rpm	
Top ceramic dimension for each position	140 × 140	180 × 180
Working time	Continue	
Positions	Two	
Stir speed	Digital screen	
Temperature range for each position	400 degree environment	
Temperature accuracy	± 1 degree	
Temperature control	PID	
Temperature display	Digital screen	
Outdoor Thermostat	To measure and control the liquid temperature value together with the front panel	
Indoor Thermostat	Controls and measures the temperature of the hot plate stirrer	
Motor power for each position	8W DC14V	8W DC14V
Product power for each position	260	410
Stirring ability for each position	50ml-1000ml	2000
Housing material	Superfine metal	
Net weight / kg)	7	10



#Letsgivemore

Thanks to your purchase, a donation will be made to a non-profit foundation that fights against breast cancer and helps most vulnerable communities.



All rights reserved © KALSTEIN France S. A. S.,
2 Rue Jean Lantier • 75001 Paris •
+33 1 78 95 87 89 / +33 6 80 76 07 10 •
<https://kalstein.eu>
KALSTEIN FRANCE, S. A. S





 **Kalstein**[®]
A different accompaniment, at your service

#LETSGIVEMORE 

**WITH THE
ACQUISITION
OF A KALSTEIN
EQUIPMENT**

YOU MAKE A CONTRIBUTION TO:

 **Fundación
Jacinto Convit**
<https://www.jacintoconvit.org/>


 **ONE
TREE
PLANTED**

<https://onetreepanted.org>

 **Humatem**
Solidarité et coopération
biomédicale
<https://www.humatem.org/>

 **MANIAPURE
FOUNDATION USA**
<https://www.maniapurefoundation.org/>

Articles Kalstein
(Social Networks and Official Websites)

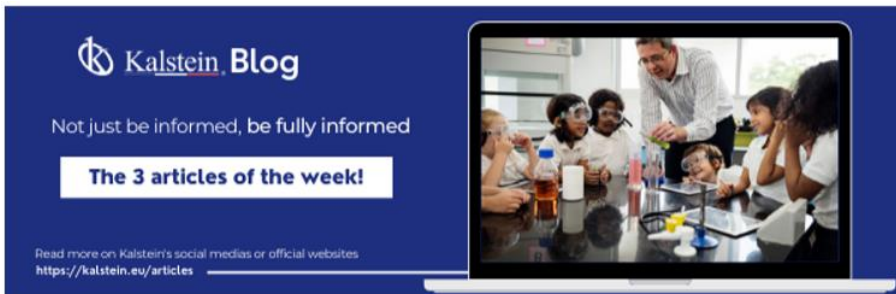

Kalstein Blog
Not just be informed, be fully informed

The 3 articles of the week!

Read more on Kalstein's social medias or official websites
<https://kalstein.eu/articles>

1
Read this ARTICLE NOW



How does the climate affect your microscope?

microscope is an instrument that allows observing and measuring small objects. This is possible through a lens system which focuses light from the object being observed into a single point.

Kalstein Blog

Article of the week:
How does the climate affect your microscope?

Read more on Kalstein's social medias or official websites
<https://kalstein.eu/articles>

2
Read this ARTICLE NOW



What is the importance of the laboratory refrigerator?

Refrigerator is a device that is one of the most common equipment that is used to maintain, in a controlled environment (refrigerated space), various fluids and substances.

Kalstein Blog

Article of the week:
What is the importance of the laboratory refrigerator?

Read more on Kalstein's social medias or official websites
<https://kalstein.eu/articles>

Kalstein Training Videos
(Youtube Channel and Official Websites)



Kalstein France - Jacinto Convit Foundation / Allianc...

23 views • 3 years ago



Kalstein's Microscopes

445 views • 1 week ago



PARIS - FRANCE
HEADQUARTER

2 Rue Jean Lantier
75001 Paris - France

Tlf: +33 (0) 01 78 95 87 89

Cel: +33 (0) 6 80 76 07 10

sales@kalstein.eu

<https://kalstein.eu/>
<https://www.kalstein.fr/>

