



LABORATORY LINES

THERMAL CYCLER



Kalstein®

A different accompaniment, at your service





Kalstein®

A different accompaniment, at your service

OUR SERVICES

Benefits and Support

In Kalstein France, we take care of the full satisfaction of our customers, that is why we provide value-added services of the highest level based on our experience.



Online Inductions and Trainings

In any part of the world, receive your induction or training from our specialized team of engineers



Quick Response

Our work team is always available to response all your consults or questions, in order to support you in any situation.



#Letsgivemore ♥

Thanks to your purchase, a donation will be made to a non-profit foundation that fights against breast cancer and helps most vulnerable communities.



Technical Support

Enjoy of personalized advice for the correct preventive and corrective maintenance of your equipment, thanks to Kalstein's manuals and articles, special catalogues and video tutorials.



Delivery Logistics

We take care of all the necessary logistics for the dispatch of your goods, whether is by sea, land or air.



Kalstein Worldwide

With more than 25 years growing with our customers, Kalstein's multiformat and modern content, is now present in more than 10 countries and increasing.





Thermal Cycler

Serie YR01867

Product Information


YR01867 PCR system allows three users to control three independent blocks in the same thermal cycler, saving your time and budget.

It combines the reliability and performance like our line of Kalstein thermal cyclers with flexible design and advanced features that adapt to your research needs.

With Android operation system, YR01867 model has been greatly optimized in terms of temperature control settings, program file management, and graphical interfaces, giving the user a more comfortable operating experience.

The model YR01867 can show its unique advantages in the case of fewer samples, more types of experiments, and more users.

The design of the independent module can help lab personnel utilize the instrument more efficiently and avoid the confusion caused by the cross usage between different users

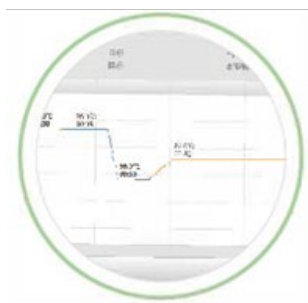

Powerful processing system:

- ✓ Advanced Android system, strong software flexibility
- ✓ 1.2GHz operating frequency, extremely fast processing speed
- ✓ 2D, 3D graphics acceleration engine, smooth running screen



Unique hot lid design:

- ✓ One-touch cover opening, convenient and quick
- ✓ Three light indicators, visible at a glance
- ✓ The spring lock design ensures that the hot lid and tube lid are in proper pressure contact



Large touch screen:

- ✓ 10.1 inch color TFT touch screen
- ✓ Familiar Android interface, comfortable software operation
- ✓ The overall status of the program is displayed intuitively



Multi-module design:

- ✓ The three modules operate independently and can be used by three people at the same time, avoiding long waits
- ✓ Parallel experiments are completed at a time, more efficient and reliable



#Letsgivemore ♥

Thanks to your purchase, a donation will be made to a non-profit foundation that fights against breast cancer and helps most vulnerable communities.



All rights reserved © KALSTEIN France S. A. S.,
Optimum Business Center 450 Rue Baden Powell,
• 34000 Montpellier, France.
Tél: +33 467158849 / +33 680760710/ +33 663810023
<https://kalstein.eu>
KALSTEIN FRANCE, S.A. S






Thermal Cycler

Serie YR01867



Three PCR programs run independently, and precise temperature control is ensured

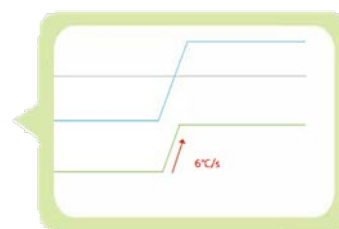

A temperature increase rate of up to 6°C/s further reduces the running time and reduces non-specific amplification.

Block forms	3×32×0.2ml	1×96×0.2ml
Temp. range	0~99.9°C	0~99.9°C
Max.heating ramp rate	6°C/sec	6°C/sec
Max.cooling ramp rate	4.5°C/sec	4.5°C/sec
Temp. uniformity	±0.2°C	±0.4°C
Temp. accuracy	±0.1°C	±0.1°C
Hot lid temp. range	20~110°C	20~110°C
Temp.adjustable rate	0.1~5°C	0.1~5°C
Gradient temp.range	1~30°C	1~30°C
Temp. control mode	Block / Calculated	Block / Calculated
Max.no. of cycle	999	999
Display	10.1 inches 1366×768 pels colored TFT touch screen	
Software	Andriod	Andriod
PC connection	Standard	Standard
Program pause function	Standard	Standard
Power down protection	Standard	Standard
Hold at 4°C	Unlimited	Unlimited
Storage	USB, SD card, On-board	USB, SD card, On-board
Communication	USB 2.0, LAN, Wifi	USB 2.0, LAN, Wifi
Power supply	AC85V-AC264V 50/60HZ	AC85V-AC264V 50/60HZ
Power input	700W	700W
Net weight	16Kg	16Kg
External dimension	502mmx394mmx293mm(L×W×H)	502mmx394mmx293mm(L×W×H)



#Let'sgivemore

Thanks to your purchase, a donation will be made to a non-profit foundation that fights against breast cancer and helps most vulnerable communities.



All rights reserved © KALSTEIN France S. A. S.,
Optimum Business Center 450 Rue Baden Powell,
• 34000 Montpellier, France.
Tit: +33 467158849 / +33 680760710/ +33 663810023
<https://kalstein.eu>
KALSTEIN FRANCE, S.A. S



Thermal Cycler

Serie YR01867



8G body storage, capacity can be extended through USB interface, and easy data import and export



Innovative design of block is easy for change. The process takes typically less than 10 seconds and requires no tools or kits.



Reserve mechanical lid opening function to ensure safe opening of hot lid in case of accident



Pull-type electromagnet open lid design, click the screen to open the lid



The temperature change rate can be set in a single section to meet demanding experimental requirements .



Simple-to-operate Android user interface is driven by the 10.1-inch colorful touch screen for better viewing and programming



Low power consumption mode and one-touch screensaver for safety and energy saving



#Letsgivemore

Thanks to your purchase, a donation will be made to a non-profit foundation that fights against breast cancer and helps most vulnerable communities.



All rights reserved © KALSTEIN France S. A. S.,
Optimum Business Center 450 Rue Baden Powell,
• 34000 Montpellier, France.
Tit: +33 467158849 / +33 680760710/ +33 663810023
<https://kalstein.eu>
KALSTEIN FRANCE, S.A. S





 **Kalstein**[®]
A different accompaniment, at your service

#LETSGIVEMORE 

**WITH THE
ACQUISITION
OF A KALSTEIN
EQUIPMENT**

YOU MAKE A CONTRIBUTION TO:



**Fundación
Jacinto Convit**

<https://www.jacintoconvit.org/>



<https://onetreepanted.org>



Humatem

Solidarité et coopération
biomédicale


<https://www.humatem.org/>



**MANIAPURE
FOUNDATION USA**

<https://www.maniapurefoundation.org/>

Articles Kalstein
(Social Networks and Official Websites)

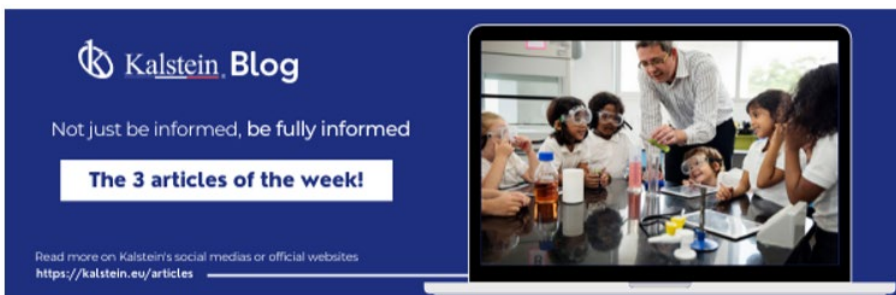

Kalstein Blog
Not just be informed, be fully informed
The 3 articles of the week!
Read more on Kalstein's social medias or official websites
<https://kalstein.eu/articles>



1
Read this ARTICLE NOW
How does the climate affect your microscope?
A microscope is an instrument that allows observing and measuring objects that are too small to be seen with the naked eye. This is achieved through a complex system of lenses and mirrors, which form...

Kalstein Blog
Article of the week:
How does the climate affect your microscope?
Read more on Kalstein's social medias or official websites
<https://kalstein.eu/articles>



2
Read this ARTICLE NOW
What is the importance of the laboratory refrigerator?
A refrigerator in a laboratory is one of the most important equipment. It is used to maintain, in a controlled environment (refrigerated space), various fluids and substances...


Kalstein Blog
Article of the week:
What is the importance of the laboratory refrigerator?
Read more on Kalstein's social medias or official websites
<https://kalstein.eu/articles>

Kalstein Training Videos
(Youtube Channel and Official Websites)



Kalstein France - Jacinto Convit Foundation / Alliances

23 views • 3 years ago

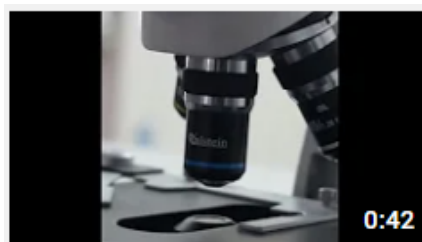

Kalstein's Microscopes

445 views • 1 week ago



MONTPELLIER - FRANCE
HEADQUARTERS

Optimum Business Center,
450 Rue Baden Powell,
34000 Montpellier - France.

Fax: +33 (0) 4 67 15 88 49

Tel: +33 (0) 6 80 76 07 10

sales@kalstein.eu

<https://kalstein.eu/>
<https://www.kalstein.fr/>

Official Websites for Latin America.

