



LABORATORY LINES

MAGNETIC STIRRERS



Kalstein®

A different accompaniment, at your service





Kalstein®

A different accompaniment, at your service

OUR SERVICES

Benefits and Support

In Kalstein France, we take care of the full satisfaction of our customers, that is why we provide value-added services of the highest level based on our experience.



Online Inductions and Trainings

In any part of the world, receive your induction or training from our specialized team of engineers



Quick Response

Our work team is always available to response all your consults or questions, in order to support you in any situation.



#Letsgivemore ♥

Thanks to your purchase, a donation will be made to a non-profit foundation that fights against breast cancer and helps most vulnerable communities.



Technical Support

Enjoy of personalized advice for the correct preventive and corrective maintenance of your equipment, thanks to Kalstein's manuals and articles, special catalogues and video tutorials.



Delivery Logistics

We take care of all the necessary logistics for the dispatch of your goods, whether is by sea, land or air.



Kalstein Worldwide

With more than 25 years growing with our customers, Kalstein's multiformat and modern content, is now present in more than 10 countries and increasing.



Thermostatic Magnetic Stirring Bath

Model YR02311/YR02312




Features

- ✓ It can be used as water bath or oil bath. Built-in strong magnetic stirrer in the bath.
- ✓ Stainless steel 304 heater.
- ✓ The magnetic stirring system drives the stirrer to rotate synchronistical, so that the solution in the bath can be heated and stirred evenly.
- ✓ DC brushless motor ensures stable performance and continuous working.
- ✓ High temperature magnet can continuously work at 300°C without losing magnetism.
- ✓ PID temperature controller ensures accurate temperature control.
- ✓ Key setting and digital display make it easy to operate.
- ✓ The model YR02311 and YR02312 are equipped with two sensors make bath temperature and container temperature display alternately.



Model YR02312

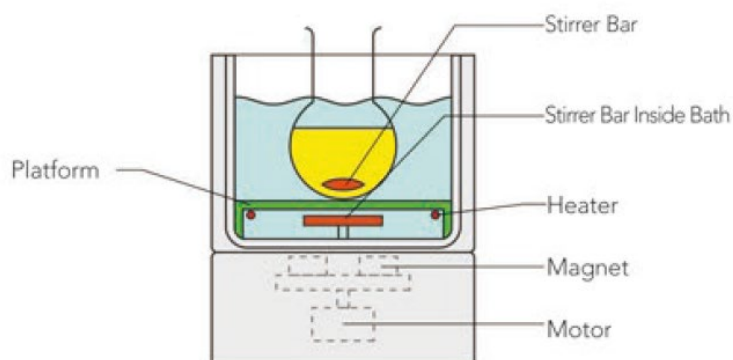

Model YR02311



#Let'sgivemore ♥

Thanks to your purchase, a donation will be made to a non-profit foundation that fights against breast cancer and helps most vulnerable communities.



All rights reserved © KALSTEIN France S. A. S.,
Optimum Business Center 450 Rue Baden Powell,
• 34000 Montpellier, France.
Tél: +33 467158849 / +33 680760710/ +33 663810023
<https://kalstein.eu>
KALSTEIN FRANCE, S.A. S



Thermostatic Magnetic Stirring Bath Model YR02311/YR02312



Model YR02311



Model YR02312

Model		
Temperature Range (°C)	Room Temperature+5 ~ 200	
Temperature Stability (°C)	±1	±1
Temperature Display	Keypad Input, Digital Display	
Stirring Speed Setting	Knob Setting	
Stirring Speed (rpm)	0~2000	0~2000
Bath Dimensions (mm)	∅ 220*110	∅ 254*130
Bath Capacity (L)	4	6.5
Heating Power (W)	500	1050
Max. Flask can be placed (ml)	3000	5000
Power Supply (V/Hz)	1~, 11 0V,60Hz or 220-240V~,50/60Hz	1~, 11 0V,60Hz or 220-240V~,50/60Hz
Dimensions (mm)	260W*280D*230H	260W*280D*260H
Net Weight (kg)	5	6



#Letsgivingmore ♥


Thanks to your purchase, a donation will be made to a non-profit foundation that fights against breast cancer and helps most vulnerable communities.



All rights reserved © KALSTEIN France S. A. S.,
Optimum Business Center 450 Rue Baden Powell,
• 34000 Montpellier, France.
Tél: +33 467158849 / +33 680760710/ +33 663810023
<https://kalstein.eu>
KALSTEIN FRANCE, S.A. S



Articles Kalstein
(Social Networks and Official Websites)

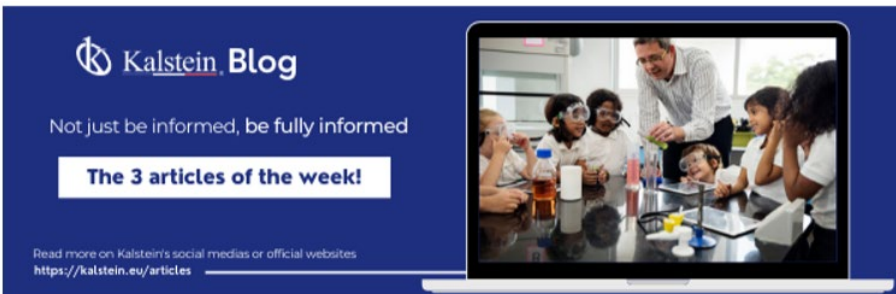

Kalstein Blog
Not just be informed, be fully informed
The 3 articles of the week!
Read more on Kalstein's social medias or official websites
<https://kalstein.eu/articles>



1
Read this ARTICLE NOW
How does the climate affect your microscope?
A microscope is an instrument that allows observing and measuring objects that are too small to be seen with the naked eye. This is achieved through a complex system of lenses and mirrors, which form...

Kalstein Blog
Article of the week:
How does the climate affect your microscope?
Read more on Kalstein's social medias or official websites
<https://kalstein.eu/articles>



2
Read this ARTICLE NOW
What is the importance of the laboratory refrigerator?
A refrigerator is a device that is one of the most common pieces of equipment used in laboratories to maintain, in a controlled environment (refrigerated space), various fluids and substances...

Kalstein Blog
Article of the week:
What is the importance of the laboratory refrigerator?
Read more on Kalstein's social medias or official websites
<https://kalstein.eu/articles>

Kalstein Training Videos
(Youtube Channel and Official Websites)



Kalstein France - Jacinto Convit Foundation / Allianc...

23 views • 3 years ago



Kalstein's Microscopes

445 views • 1 week ago



MONTPELLIER - FRANCE
HEADQUARTER
Optimum Business Center,
450 Rue Baden Powell,
34000 Montpellier - France.

Fax: +33 (0) 4 67 15 88 49
Tif: +33 (0) 6 80 76 07 10
sales@kalstein.eu

<https://kalstein.eu/>
<https://www.kalstein.fr/>

Official Websites for Latin America.

