




LABORATORY LINES

VACUUM OVEN



Kalstein®

A different accompaniment, at your service

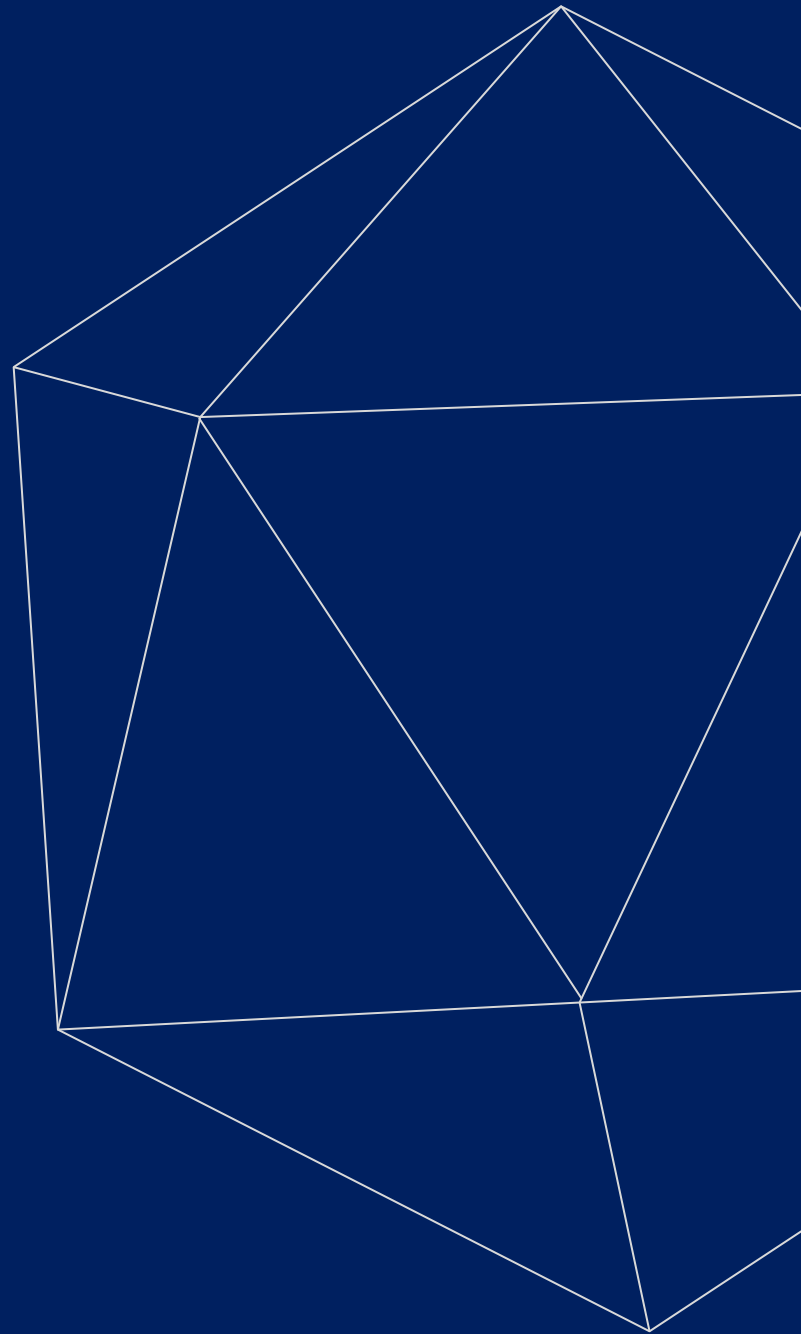



Kalstein®

A different accompaniment, at your service

OUR SERVICES

Benefits and Support

In Kalstein France, we take care of the full satisfaction of our customers, that is why we provide value-added services of the highest level based on our experience.



Online Inductions and Trainings

In any part of the world, receive your induction or training from our specialized team of engineers



Quick Response

Our work team is always available to response all your consults or questions, in order to support you in any situation.



#Letsgivemore ❤️

Thanks to your purchase, a donation will be made to a non-profit foundation that fights against breast cancer and helps most vulnerable communities.



Technical Support

Enjoy of personalized advice for the correct preventive and corrective maintenance of your equipment, thanks to Kalstein's manuals and articles, special catalogues and video tutorials.



Delivery Logistics

We take care of all the necessary logistics for the dispatch of your goods, whether is by sea, land or air.



Kalstein Worldwide

With more than 25 years growing with our customers, Kalstein's multiformat and modern content, is now present in more than 10 countries and increasing.



Vacuum drying oven Serie YR



Kalstein's Vacuum Drying oven widely used in biochemical, chemical and pharmaceutical, health care, agricultural research, environmental protection and other research applications for powder Drying and roasting, as well as disinfection and sterilization of glass containers used.

Technology advantage:

- ✓ The vacuum environment greatly reduces the boiling point of the liquid to be removed, so vacuum drying can be easily applied to heat sensitive substances. For samples that are not easy to dry, such as powder or other granular samples, vacuum drying can effectively shorten the drying time.

Structure features

- ✓ The inner chamber is of arc-shaped design, the outer is made of cold-rolled steel plate, and the surface is electrostatically sprayed,
- ✓ Stainless steel interior or cold-rolled steel interior.
- ✓ The tightness of the cabinet door can be adjusted by the customer. The integrally formed silicone rubber door sealing ring ensures the high vacuum in the cabinet
- ✓ The control panel adopts microcomputer intelligent PID control, with timing and over temperature alarm functions.
- ✓ Tempered and explosion-proof double-layer glass door.



YR05988 (S) -YR05991 (S) Serie: Individual shelf heating



YR05262 (SS) - YR05265-1 (S) series: 2 sided heating
YR05992 (S) - YR05995 (S) series: 4 sided heating



#Letsgivemore ♥

Thanks to your purchase, a donation will be made to a non-profit foundation that fights against breast cancer and helps most vulnerable communities.



All rights reserved © KALSTEIN France S. A. S.,
2 Rue Jean Lantier • 75001 Paris •
+33 1 78 95 87 89 / +33 6 80 76 07 10 •
<https://kalstein.eu>
KALSTEIN FRANCE, S. A. S



Vacuum drying oven
Serie YR



#Letsgivemore

Thanks to your purchase, a donation will be made to a non-profit foundation that fights against breast cancer and helps most vulnerable communities.





All rights reserved © KALSTEIN France S. A. S.,
2 Rue Jean Lantier • 75001 Paris •
+33 1 78 95 87 89 / +33 6 80 76 07 10 •
<https://kalstein.eu>
KALSTEIN FRANCE, S. A. S



Vacuum drying oven Serie YR





DZF-60 series heating diagram

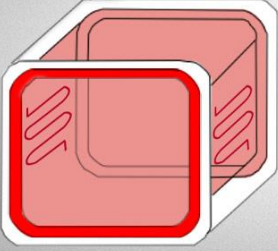

- Two sides heating technology(DZF-60)
- Stainless steel inner (type with "s")

DZF-50 series heating diagram




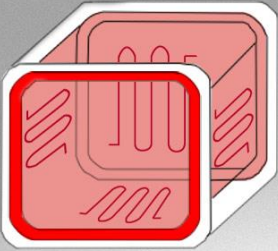

- Four sides heating technology(DZF-50)
- Stainless steel inner,qualified aluminum sheet

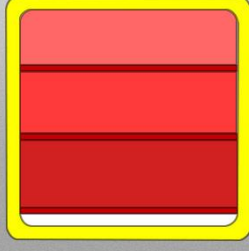
DZF-30 series heating diagram


- Independent temperature control technology(DZF-30)
- Stainless steel inner,qualified aluminum sheet

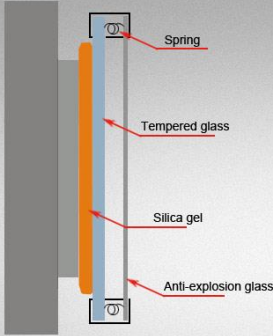
- Anti pressure protection device
- LED lamp




1.Vacuum pump

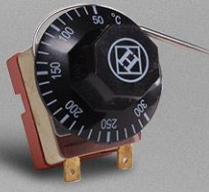


2.Independent temperature controller



3.Tray material, LED lighting lamp



4. 30 sections temperature controller

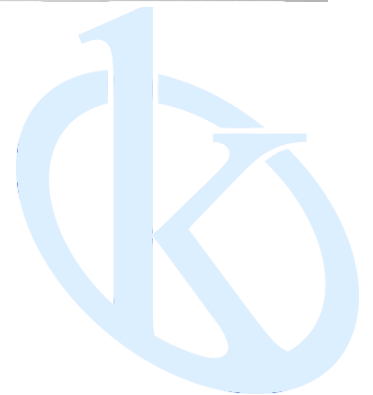

#Letsgivemore 

Thanks to your purchase, a donation will be made to a non-profit foundation that fights against breast cancer and helps most vulnerable communities.



All rights reserved © KALSTEIN France S. A. S.,
2 Rue Jean Lantier • 75001 Paris •
+33 1 78 95 87 89 / +33 6 80 76 07 10 •
<https://kalstein.eu>
KALSTEIN FRANCE, S. A. S



Vacuum drying oven Serie YR



Technical Details

Model	Voltage V	Temp. range°C	Timer	Accuracy °C	Vacuum degree Pa	Power KW	Interior (mm) (H×W×D)	Exterior (mm) (H×W×D)	Overall (mm) (H×W×D)	N.W KG
YR05262 (SS)	220	50~200	0~999	±1	< 133	0.8	280×300×300	450×635×465	570×720×545	42.5
YR05263 (SS)						1.2	400×415×350	570×750×515	690×830×600	65
YR05265 (SS)						2	485×440×425	620×780×600	740*900*710	80
YR05265-1 (S)						2	485×440×500	1580×700×720	1770×790×790	120
YR05988 (S)						0.8	315×265×280	450×635×465	570×720×545	45
YR05989 (S)						1.2	400×370×350	570×750×515	690×830×600	68
YR05990 (S)						2	485×400×425	620×780×600	740*900*710	85
YR05991						2	485×440×500	1580×700×720	1770×790×790	120
YR05992 (S)						1	280×265×245	450×635×465	570×720×545	43
YR05993 (S)						1.4	370×370×315	570×750×515	690×830×600	66
YR05994 (S)						2	455×440×390	620×780×600	740*900*710	81
YR05995 (S)						2	485×440×500	1580×700×720	1770×790×790	120

The modes with "S" is stainless steel inner

Optional parts

- 1、 vacuum pump
- 2、 individual temperature limit



#Letsgivemore

Thanks to your purchase, a donation will be made to a non-profit foundation that fights against breast cancer and helps most vulnerable communities.



All rights reserved © KALSTEIN France S. A. S.,
2 Rue Jean Lantier • 75001 Paris •
+33 1 78 95 87 89 / +33 6 80 76 07 10 •
<https://kalstein.eu>
KALSTEIN FRANCE, S. A. S





#LETSGIVEMORE

**WITH THE
ACQUISITION
OF A KALSTEIN
EQUIPMENT**

YOU MAKE A CONTRIBUTION TO:



Fundación
Jacinto Convit

<https://www.jacintoconvit.org/>



<https://onetreepanted.org>



Humatem

Solidarité et coopération
biomédicale

<https://www.humatem.org/>



**MANIAPURE
FOUNDATION USA**

<https://www.maniapurefoundation.org/>

Articles Kalstein
(Social Networks and Official Websites)



Kalstein Blog
Not just be informed, be fully informed

The 3 articles of the week!

Read more on Kalstein's social medias or official websites
<https://kalstein.eu/articles>

1
Read this ARTICLE NOW



How does the climate affect your microscope?

microscope is an instrument that allows observing and measuring small objects. This is possible through a series of lenses that focus light on the object to be observed.

Kalstein Blog

Article of the week:
How does the climate affect your microscope?

Read more on Kalstein's social medias or official websites
<https://kalstein.eu/articles>

2
Read this ARTICLE NOW



What is the importance of the laboratory refrigerator?

Refrigerator is a device that is one of the most common equipment that is used to maintain, in a controlled environment (refrigerated space), various fluids and substances.

Kalstein Blog

Article of the week:
What is the importance of the laboratory refrigerator?

Read more on Kalstein's social medias or official websites
<https://kalstein.eu/articles>

Kalstein Training Videos
(Youtube Channel and Official Websites)



Kalstein France - Jacinto Convit Foundation / Allianc...
23 views • 3 years ago



Kalstein's Microscopes
445 views • 1 week ago



PARIS - FRANCE
HEADQUARTER
2 Rue Jean Lantier
75001 Paris - France

Tlf: +33 (0) 01 78 95 87 89
Cel: +33 (0) 6 80 76 07 10
sales@kalstein.eu

<https://kalstein.eu/>
<https://www.kalstein.fr/>

