




# LABORATORY LINES

LAMINAR FLOW CABINET



**Kalstein**®

A different accompaniment, at your service

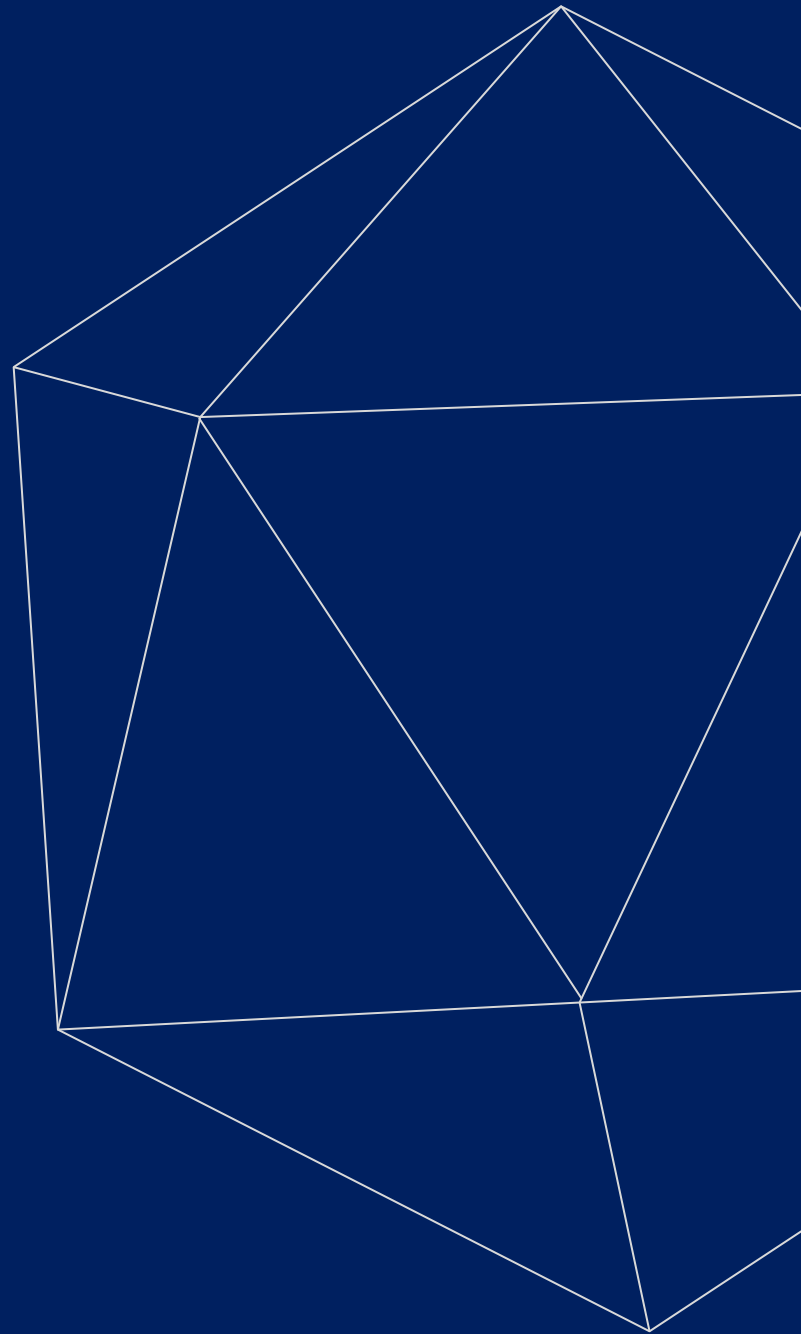



# Kalstein®

A different accompaniment, at your service

OUR SERVICES

## Benefits and Support

In Kalstein France, we take care of the full satisfaction of our customers, that is why we provide value-added services of the highest level based on our experience.



### Online Inductions and Trainings

In any part of the world, receive your induction or training from our specialized team of engineers



### Quick Response

Our work team is always available to response all your consults or questions, in order to support you in any situation.



### #Letsgivemore ♥

Thanks to your purchase, a donation will be made to a non-profit foundation that fights against breast cancer and helps most vulnerable communities.



### Technical Support

Enjoy of personalized advice for the correct preventive and corrective maintenance of your equipment, thanks to Kalstein's manuals and articles, special catalogues and video tutorials.



### Delivery Logistics

We take care of all the necessary logistics for the dispatch of your goods, whether is by sea, land or air.



### Kalstein Worldwide

With more than 25 years growing with our customers, Kalstein's multiformat and modern content, is now present in more than 10 countries and increasing.



## KALSTEIN® Laminar Flow Cabinet

### Laminar Flow Cabinet (Double person)

#### Model Series YR05267-1



#### Product Description

YR series (double people) laminar flow cabinet combined with perfect airflow mechanical design for the operator to create a partial purification work environment so that experiments can be carried out smoothly.

Applicable in the field of protection of biotechnology, pharmaceutical, life sciences and precision electronic products..

#### Product Introduction

- Integrated LCD dot matrix LCD display, real-time display the status of lighting, sterile, air flow and other parameters.
- User-friendly design, with sterilization, lighting, interlock functions to protect the operator from UV light burns and ensure safe operation.
- Wind speed multi-level adjustable, the air volume can be adjusted according to the user.

#### Standard Features

- ✓ CLEANFLOW™ clean air technology
- ✓ Centrifugal wind wheel low-power motor, based on airflow mechanics principle, arc-shaped large-diameter inlet design to reduce turbulence generated to ensure low noise operation of the system.
- ✓ Reasonable duct design, optimized net press box, with micro-mesh design layer outflow wind deflector, convex, no dead ends, to ensure the uniformity and stability of the airflow.
- ✓ Standard early-efficient and efficient double-filter, primary filter to avoid airflow dust particles to enter and extend the life of high efficiency filter. High efficiency filter without clapboard HEPA, retention efficiency of 99.99%.
- ✓ According to the customer's requirements, CLEANFLOWseries with the vertical flow and horizontal airflow.



#Letsgivemore ❤️

Thanks to your purchase, a donation will be made to a non-profit foundation that fights against breast cancer and helps most vulnerable communities.



All rights reserved © KALSTEIN France S. A. S.,  
2 Rue Jean Lantier • 75001 Paris •  
+33 1 78 95 87 89 / +33 6 80 76 07 10 •  
<https://kalstein.eu>  
KALSTEIN FRANCE, S. A. S



## KALSTEIN® Laminar Flow Cabinet

### Laminar Flow Cabinet (Double person)

#### Model Series YR05267-1



#### Technical Parameters

Model	YR05267-1
Cleanliness Rating	The international ISO 14644.1 standard, CLASS 5 (U.S. Federal 209E, 100 level)
Wind Speed	0.3~0.6m/s(Adjustable)
Noise (dB)	≤62
Vibration half peak	≤3μm (X、Y、Z direction)
High efficiency filter	Efficiency 99.99%
Illumination (LX)	≥300
Work nvironment	Ambient temperature 10~30°C, Relative humidity below 70%, No obvious vibration and dust around
Overall size (mm)(H×W×D)	1585×1460×735
Inner size (mm)(H×W×D)	500×1300×550
Voltage(V)	AC220V/50Hz
Power (W)	400
Net weight (KG)	210
Shipping weight(KG)	235
Packing size(mm) (H×W×D)	1030×1650×855



#### #Let'sgivemore ♥

Thanks to your purchase, a donation will be made to a non-profit foundation that fights against breast cancer and helps most vulnerable communities.



All rights reserved © KALSTEIN France S. A. S.,  
2 Rue Jean Lantier • 75001 Paris •  
+33 1 78 95 87 89 / +33 6 80 76 07 10 •  
<https://kalstein.eu>  
KALSTEIN FRANCE, S. A. S







 **Kalstein**<sup>®</sup>  
A different accompaniment, at your service

#LETSGIVEMORE 

**WITH THE  
ACQUISITION  
OF A KALSTEIN  
EQUIPMENT**

YOU MAKE A CONTRIBUTION TO:



Fundación  
Jacinto Convit

<https://www.jacintoconvit.org/>



<https://onetreepanted.org>



**Humatem**

Solidarité et coopération  
biomédicale

<https://www.humatem.org/>



**MANIAPURE  
FOUNDATION USA**

<https://www.maniapurefoundation.org/>



**Articles Kalstein**  
(Social Networks and Official Websites)



**Kalstein Blog**  
Not just be informed, be fully informed

**The 3 articles of the week!**

Read more on Kalstein's social medias or official websites  
<https://kalstein.eu/articles>



**1**  
Read this ARTICLE NOW

**How does the climate affect your microscope?**

**Kalstein Blog**

**Article of the week:**  
How does the climate affect your microscope?

Read more on Kalstein's social medias or official websites  
<https://kalstein.eu/articles>



**2**  
Read this ARTICLE NOW

**What is the importance of the laboratory refrigerator?**

**Kalstein Blog**

**Article of the week:**  
What is the importance of the laboratory refrigerator?

Read more on Kalstein's social medias or official websites  
<https://kalstein.eu/articles>

**Kalstein Training Videos**  
(Youtube Channel and Official Websites)



**Kalstein France - Jacinto Convit Foundation / Allianc...**  
23 views • 3 years ago



**Kalstein's Microscopes**  
445 views • 1 week ago



PARIS - FRANCE  
**HEADQUARTER**  
2 Rue Jean Lantier  
75001 Paris - France

**Tlf:** +33 (0) 01 78 95 87 89  
**Cel:** +33 (0) 6 80 76 07 10  
[sales@kalstein.eu](mailto:sales@kalstein.eu)

<https://kalstein.eu/>  
<https://www.kalstein.fr/>

